



How to Disconnect a Philips Respironics V60

This guide will show you how to disconnect the power of a Philips Respironics V60 ventilator.

Written By: Rie Kuraoka

INTRODUCTION


Use this guide to properly disconnect the power of the Philips Respironics V60 ventilator from a power source.

Step 1 — shut down the ventilator



- Press and release the **ON/Shutdown** key, opening the **Shutdown** window in the screen.
 - Select **Ventilator Shutdown**, which will shut down the ventilator.
- ⚠ Improper shutdown may cause a **Power has been restored** message the next time the ventilator is turned on.

Step 2 — disconnect from power

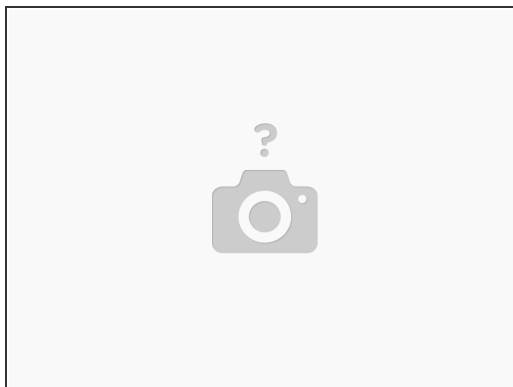


- Disconnect the ventilator from the AC power.

Step 3 — How to Disconnect a Philips Respironics V60



- Unplug the ventilator

Step 4

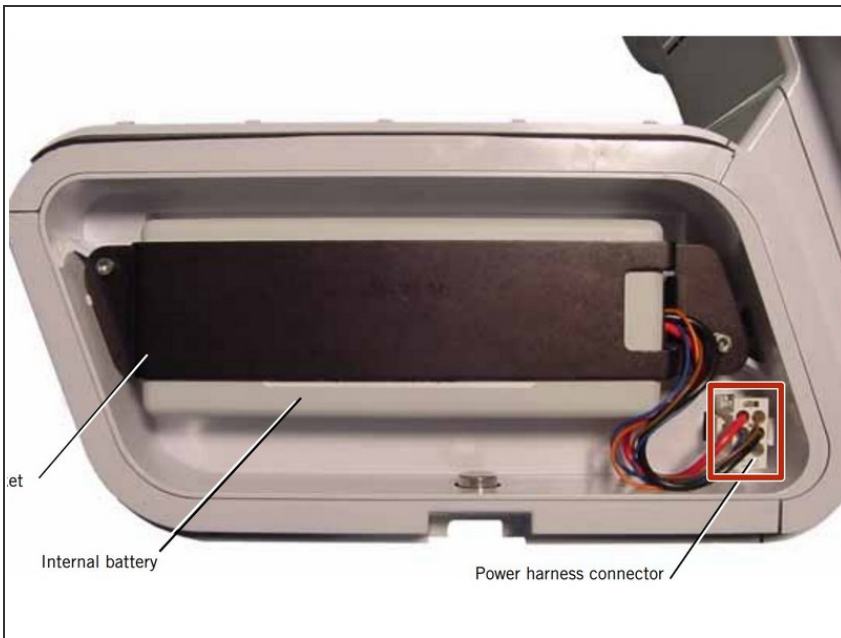

- Disconnect the internal battery from the power harness connector.

To reassemble your device, follow these instructions in reverse order.