



# Hamilton G5 Rail Removal

This guide will take you through the process of how to remove the rail on a Hamilton G5.

Written By: Sierra Saenz



---

## INTRODUCTION


For more information on removing the rail, please refer to **section 10** in the service manual.

---

### TOOLS:


- [Allen Wrench](#) (1)
-

## Step 1 — Rail

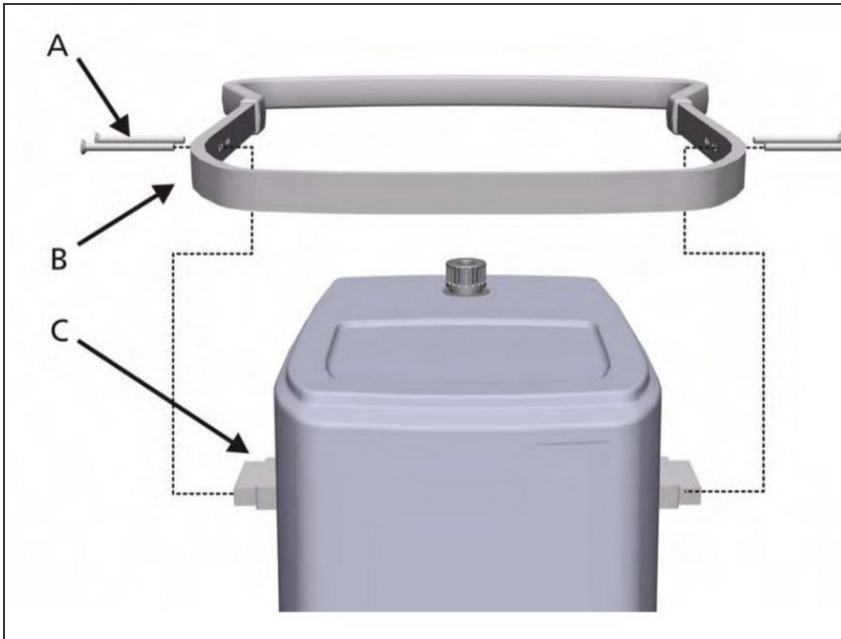


- Remove 4 Hex Allen screws (A) from the Rail surrounding the Ventilation Unit.

## Step 2

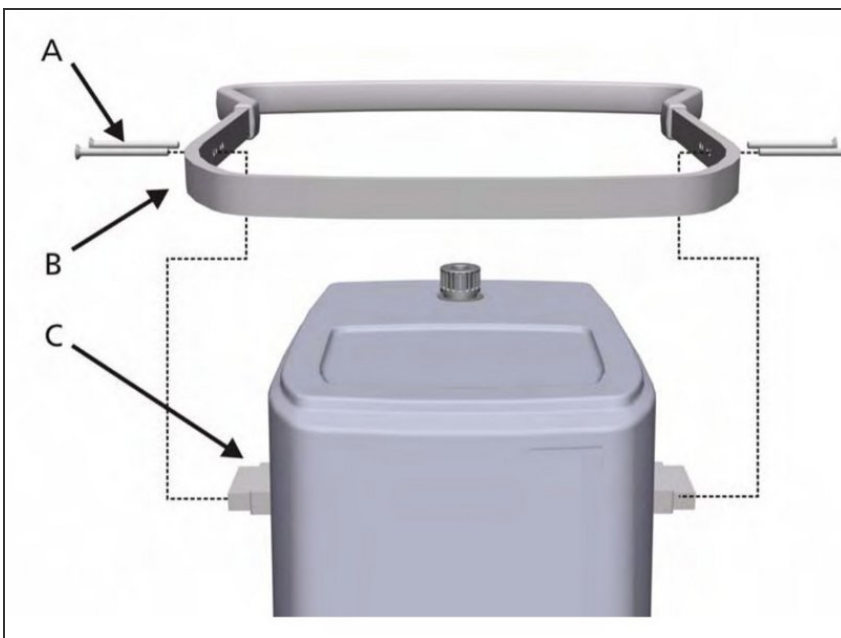


- Remove the Rail (B).

### Step 3

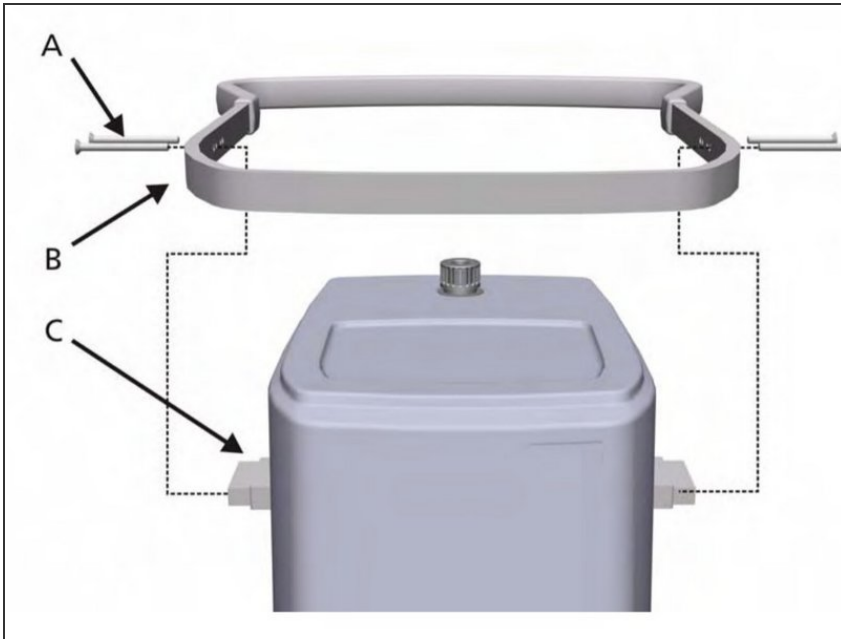

- The Standoff Brackets (C) remain with the Main Cover.

### Step 4

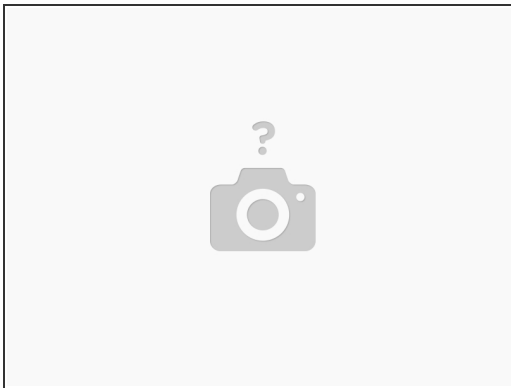

- Place aside in a safe location.

To reassemble your device, follow these instructions in reverse order.